



GLOBAL HEALTHCARE

Superior supplies for exceptional patient care

Medical Inkjet Printer For X Ray Film JD-P5670

Feature and usage:

Medical inkjet printer is an economical image printer developed for the needs of the medical image output field. This printer was composing by printer and CISS. Cooperate with special medical image ink to achieve high-precision, high-resolution and high-gray printing output of medical images. For the big CISS ink tank, one time add ink, could print max 2000 sheets 8x10 size film.

> **Professional:** For medical imaging characteristics, tailor-made, patented technical quality, beautiful appearance, delicate touch, high resolution, intuitive diagnosis.

> **Efficacy:** Avoid adding ink frequently to improve work efficiency. The large-capacity ink cartridge adopts a large opening design, which can easily add ink.

Environmental protection: PET material, environmental protection concept, green products.

Parameter: Print image principle	Inkjet
Resolution	9600x2400
Print speed	8x10 60sheets/h, 13x17 20 sheets/h
Tray volume	Max 100 sheets
Ink tank volume	500/1000 ml
Color Support	Black & White, Color image output
Ink color support	5 color
Compatible film	Blue color film, White color film
Weight	17KG
Operate temperature	0-35℃
Operate humidity	35%-75%





GLOBAL HEALTHCARE

Superior supplies for exceptional patient care



- **Usage**

Medical ink was special used for medical image output. Epically used for CR, DR, CT, MRI image output.

- **Feature**

Original ink:

The print image seems reddish when use the *original ink* that provide by Canon, Epson company.

Medical ink:

Medical ink was all black color which with different black density. To form the gray scale contrast.

Ink details for different type printer

Printer Brand	Color Support	Detail	Recommend print model
Canon G2000, etc.	4	C M Y K	Canon G1010,
Canon	5	C M Y K BK	Ix6580, 6780, etc.
Epson	4	C M Y K	L310, L1300, etc.
Epson	6	C M Y K LC LM	L805, L1800, etc.
Note: Similar model also available for print film.			

Packing details:

Storage

Please store the film in cool and dry place without dust. Keep away from direct sunlight.

Packing detail


Printer	Brand	Package
Canon		500 ml/bottle
Epson		70 ml/bottle 500 ml/bottle



GLOBAL HEALTHCARE

Superior supplies for exceptional patient care

Print contrast with different ink on blue film & white film

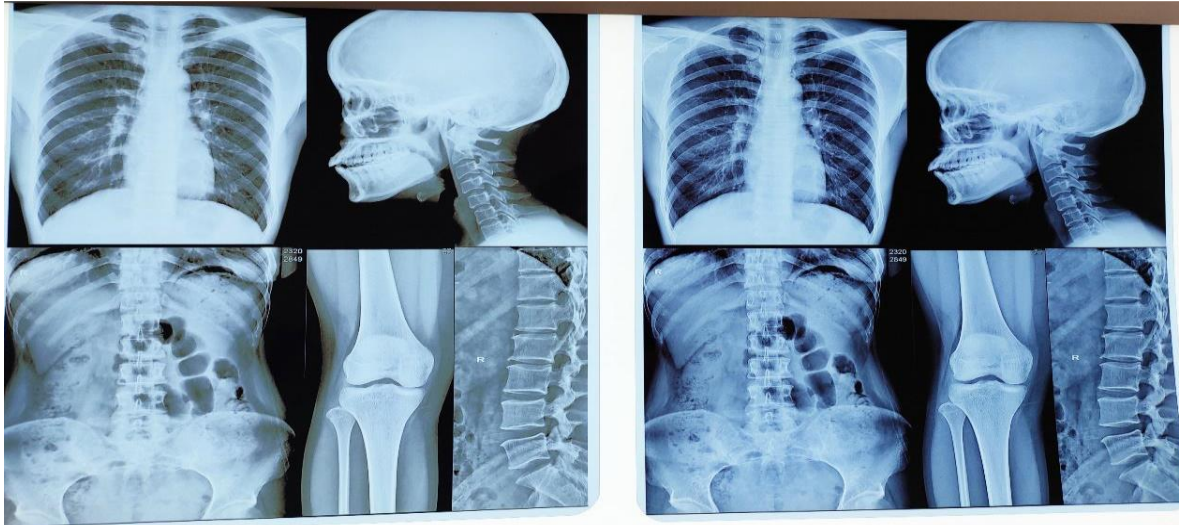

Blue base film: **Left**, Original ink print effect. **Right**, Medical ink print effect.