




Cordis 3

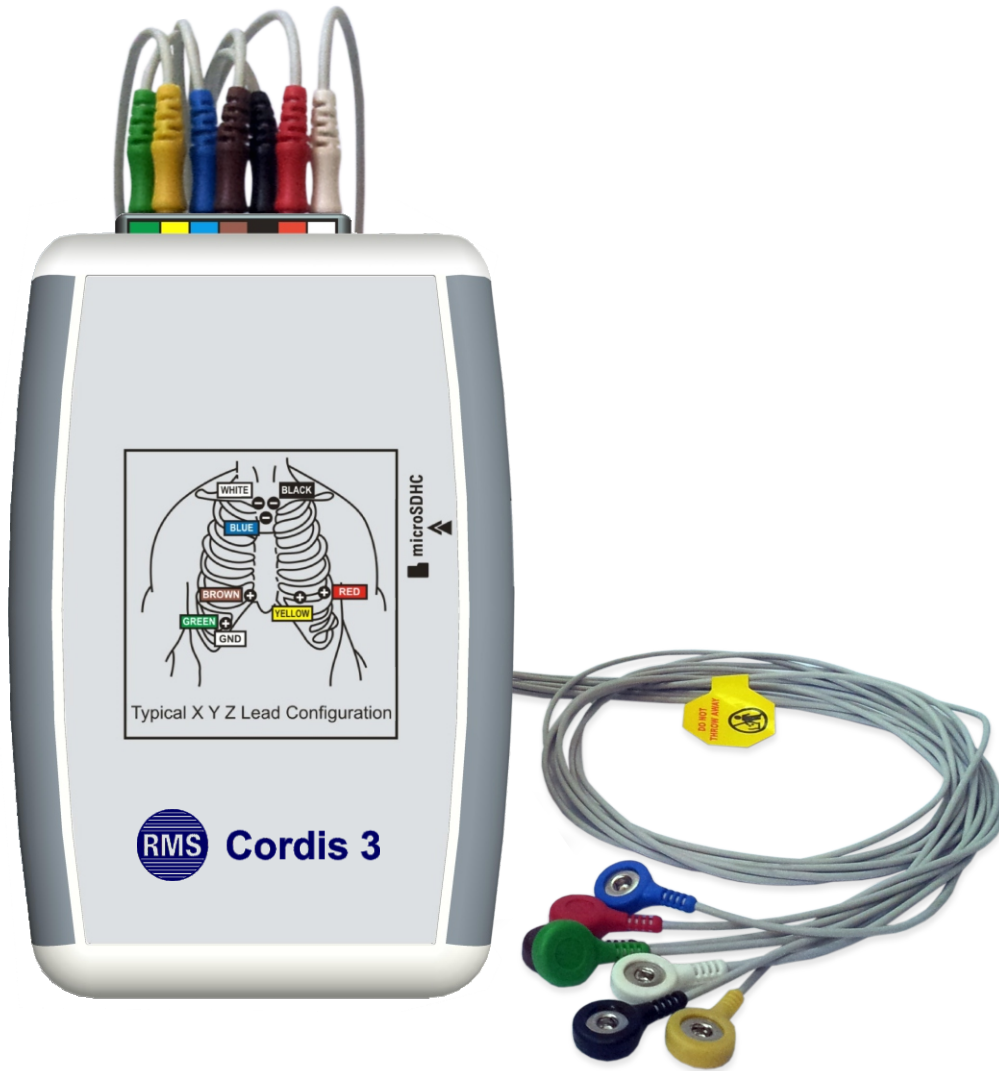

3 CHANNEL ECG HOLTER

Digital Holter Recorder

Digital Holter Recorder is a three-channel, seven-lead solid-state device designed for 24-hour continuous recording of ambulatory ECG data. Its nonvolatile, removable, digital storage and rugged, compact design ensure reliable operation and long life, while its light weight enables patients comfort during daily activities.

Cordis 3 - Technical Specifications

RECORDER

Recording Channels	3 Channels
Leads	7
Storage	SDHC card (4GB)
Storage Indication	LED indication
Event Marking	Event button
Recording	Full information without compressing
Duration	24hrs
Dimensions	120mm X 77mm X 24mm
Power Supply	Alkaline battery 1.5V (AA) LR6 - 2 nos.

ANALYSER

- 3 lead arrhythmia analysis
- Multiple Template editing to analyze VE, SVE, etc.
- Heart rate variability with 3D Spectral
- Hourly count menu
- Atria fibrillation / Flutter Scanning
- Pacemaker auto detection and analysis
- Search function based on arrhythmia type
- QT analysis
- ST analysis with Reset
- ASCII Output

MINIMUM COMPUTER CONFIGURATION

OS: Windows XP Pro/Windows 7 Professional 32 bit, Processor: Core2Duo or higher, RAM: 2GB or higher, 40 GB harddisk or higher, CD/DVD Optical Drive, Screen Resolution 1024x768 or higher




RECORDER POUCH



**SHOULDER BELT
WAIST BELT**



ELECTRODE

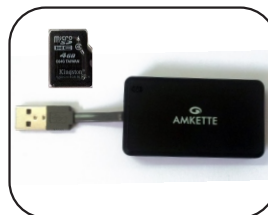

**SDHC CARD
AND READER**



PATIENT CABLE

Remote Service Program

Online technical support is available for RMS Head office. Customer has to provide internet connection.

Since R&D is a continuous process, to make changes in product features, specifications, aesthetics and/or to discontinue the same at any time without notice or obligation

Certified ISO 9001:2015, EN ISO 13485:2016 company

Recorders & Medicare Systems (P) Ltd.

Plot No. 11, Information Technology Park, HSIIDC, Sector 22, Panchkula-134114,
Haryana, INDIA, Ph.: +91-172-2564196, 2565196, Toll Free: 1800-1200-767,
www.rmsindia.com, sales@rmsindia.com
Unit II: Plot No. 9, Industrial Area, Phase 1, Panchkula-134113, Haryana, INDIA

